
nuvo

Pure 1

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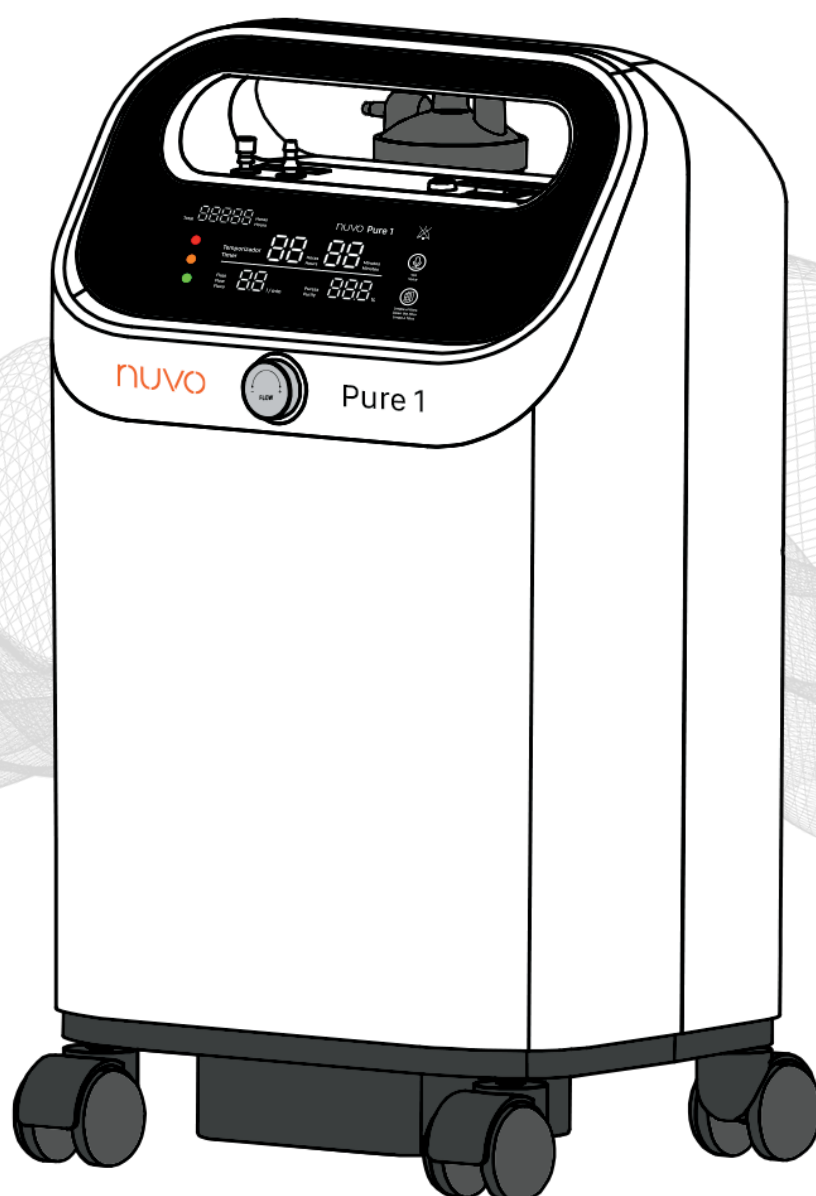1. Safety notes

This equipment is NOT a life-support device.

If the patient is unable to communicate discomfort while using the equipment, additional supervision is required.

Glossary of Symbols

	ON Power on		OFF Power off
	Prohibited action		Mandatory action
	Warning		Consult the instructions for use
	Manufacturer Name and address		Date of manufacture
	Serial number		Unique Device Identifier (UDI)
	Alternating current (AC)		Class II equipment
	Applied part type BF		Indicates medical device
	No smoking		No open flames Fire/ignition sources
	Do not sit		Do not push
	Air outlet		Non-ionizing electromagnetic radiation
	Cancel/mute alarm		Temperature limit
	Humidity limit		Atmospheric pressure limit
	Keep away from rain		Fragile Handle with care
	This side up		Stacking limit by number
	Recyclable/recovery		Do not dispose as household waste Send for separate collection and recycling
	Authorized representative in the European Community		IP21: Protected against contact with fingers or objects >12.5 mm and against vertical dripping water (no harmful effect)
	This product has been evaluated by the manufacturer and complies with EU safety, health, and environmental protection requirements.		

Note (electrical hazard): Incorrect operation may cause electric shock. **DO NOT disassemble the equipment.** Service must be performed by qualified personnel only.

1.2 Before installation

Operating environmental conditions

- **Ambient temperature:** +5 °C to +40 °C
- **Relative humidity:** 20% to 80% (non-condensing)
- **Atmospheric pressure:** 86 kPa to 106 kPa
- The device can operate up to 2000 m altitude at 21 °C without performance degradation.
- Performance may be affected outside these ranges.

Note: if the storage/transport temperature was **below 5 °C**, allow the device to remain at normal operating temperature for **more than 4 hours before use**.

Location and environment

DO NOT use in explosive atmospheres.

To reduce fire/explosion risk, keep away from heat sources, ignition sources, solvents, aerosols, etc.

- Place the device in a well-ventilated area, free from contaminants/smoke, protected from weather, and with adequate lighting.
- Position the device so that the power cord and tubing/cannula do not create a tripping hazard; the power cord must remain accessible for disconnection.

Electrical power

Use only the supplied power cord.

Verify that the power outlet specifications match the device rating label (rear panel).

Electrical versions: AC 110 V~ 50–60 Hz or AC 220 V~ 50–60 Hz (depending on version; check the rating label).

If mains voltage is unstable or exceeds the version range, install a voltage stabilizer before use.

The device requires electrical power to operate. For continuity during power outages, a backup power source is recommended.

DO NOT use extension cords or adapters (risk of sparks and fire).

Other precautions

DO NOT use in magnetic environments (MRI, X-ray, etc.): may cause malfunction.

DO NOT block filters or ventilation openings with sheets, curtains, or other objects.

Keep pets, pests, and children away from the device; prevent inhalation or ingestion of small objects.

Risk of strangulation from cables and hoses, especially if excessively long.

Prohibitions

DO NOT modify the equipment.

DO NOT connect accessories or devices not specified in this manual.

1.3 Fire or explosion warning

-  Petroleum-based lubricants, oils, lotions, and cosmetics may be flammable; their use while the device is operating is dangerous.
-  This device delivers high-concentration oxygen, which accelerates combustion. Keep tobacco, cigarettes, and open flames away from the device.

1.4 Measures to reduce risks

-  **Do not** place or store the device where water or other liquids may drip onto it.
If the device falls into water: **DO NOT touch it**, immediately disconnect the power supply, and contact your supplier.
-  Once connected to the power supply, do not leave the device unattended.
If you feel the oxygen flow is insufficient, contact your supplier or physician and adjust the flow according to their instructions.
Use of non-recommended accessories or spare parts may reduce performance and void manufacturer liability.
If discomfort or a medical emergency occurs during therapy, seek medical attention immediately.
Avoid any sparks near the device (including static electricity from friction).
If the device malfunctions or is damaged (device, cable, or plug), contact the supplier for inspection and repair.
-  **Do not** move the device without first disconnecting the power supply.
Do not drip or insert substances into device openings.
Service and maintenance are prohibited while the medical device is in use.
-  Elderly patients or others unable to express discomfort may require additional supervision to prevent injury.

1.5 Precautions

-  The patient may be the intended operator. Operating, care, and safety instructions apply to the operator.
Sale restricted by or on the order of a physician.
Intended for trained operators or those experienced with similar equipment.
Cleaning and disinfection may be performed by the patient.
-  **Do not** use for more than 24 continuous hours and not suitable for critically ill patients.
Report serious incidents related to the device to the manufacturer and the competent authority in your country.
-  **Contraindications:** not suitable for surgery, emergencies, or critically ill patients. Contraindicated for patients with oxygen toxicity or oxygen allergy.

Possible side effects during treatment:

- Dry mouth, nose, and throat
- Dry or cold nose
- Runny nose
- Nasal congestion or cold symptoms
- Itching in the nasal cavity

1.6 Service life

- **Oxygen concentrator:** 5 years
- **Nasal cannula (single-use, sterilized with ethylene oxide):** expiration and production dates indicated on packaging

1.7 Operating principle

- 1 The device draws air through the external intake filter. The filtered air enters the compressor (via a suction resonator and fine filter). The pressurized air passes through a heat exchanger that reduces its temperature. A valve system directs the air to one of two molecular sieve beds.
- 2 The sieve adsorbs nitrogen from the air, allowing oxygen-enriched gas to pass to the pressure regulator. While one bed produces product gas, the other is purged to release the adsorbed nitrogen; this process is called Pressure Swing Adsorption (PSA).
- 3 The flow delivered to the patient is then adjusted using the flow meter/valve.
- 4 Finally, the gas passes through a fine filter and an oxygen concentration sensor before exiting through a fire-resistant outlet.

1.8 Basic information

Intended use

This oxygen concentrator is used to continuously provide oxygen-enriched product gas with a purity of 90%–96% (V/V) for patients with health conditions that cause low blood oxygen levels (hypoxemia). It is intended for adult patients in home or hospital environments.

Attention: this is a medical device. Read the instructions carefully before operation. If in doubt, consult a healthcare professional.

Intended user (operator)

- **Age:** Adult (over 18 years).
- **Minimum knowledge:** Ability to read and understand text and numbers; read this manual.
- **Language:** Spanish, Portuguese or English
- **Experience:** Training provided via video.
- **Allowed limitations:** Mild visual impairment for reading or corrected visual acuity up to logMAR 0.5 (6/19 or 20/63) or better.
- **Hearing impairment:** Up to 40% impairment resulting in 60% normal hearing between 50 Hz and 2 kHz.

Patients and caregivers must meet equivalent requirements.

Intended environment

- Home and healthcare environments.

Intended patient population

- Patients requiring supplemental oxygen.

Indications

- Suitable for adult patients who require supplemental oxygen to treat or relieve their condition.

Contraindications: this device must not be used as a life-support system.

Clinical benefits

- Home oxygen concentrators provide continuous and stable oxygen therapy support to help patients with chronic diseases improve physiological indicators, slow disease progression, and improve quality of life.

2. Product features

Product structure

This product consists of:

- Oxygen concentrator · 1 unit
- Nebulization kit · 1 unit
- Nasal cannula · 1 unit
- Firebreak valve · 1 unit
- Humidifier bottle · 1 unit
- Silicone connection tube · 1 unit
- First-stage black filter (pre-filter) · 2 units
- Second-stage white HEPA filter · 1 unit
- Remote control · 1 unit


Structural features

- Full ABS plastic housing.
- Accumulated time function (total operating time displayed on LED screen).
- Scheduled shutdown (timer) function.
- Compressor includes a safety valve.
- Compressor includes thermal protection.
- Power failure alarm function.

Main technical parameters

Oxygen performance and output

Item	Specification	Notes
Recommended maximum flow	5 L/min	
Flow range	0,5 L/min~5 L/min	When the nominal outlet pressure is 7 kPa
Flow variation	0,5 L/min	At recommended max flow and 7 kPa back pressure
Oxygen purity (O ₂)	90%~96% (V/V)	When nominal outlet pressure is zero (fig. 1)
Outlet pressure	0,02 MPa ~ 0,07 MPa	
Oxygen outlet temperature	< 46 °C	
Minimum operating time	≥ 30 min	
Maximum device output	7 L/min	Above 5 L/min oxygen purity drops below 90%

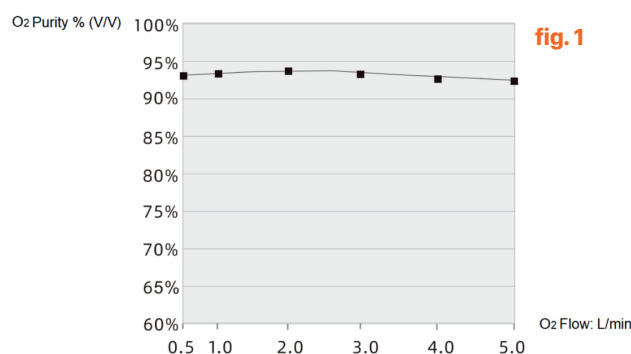! Refer to the technical label on the rear of the device to verify the actual electrical supply.

Safety, Alarms, and Functions

Item	Specification	Notes
Safety valve	250 kPa $\pm$ 30 kPa	At compressor outlet
Low purity	Visual and audible alarm; continues operating	
Low voltage	Visual and audible alarm; continues operating	Threshold depends on electrical version
Power failure	Audible alarm; stops operation	
Noise	$\leq$ 48 dB(A)	

Power supply and consumption


Item	Specification	Notes
Power supply	CA 110 V~ 50–60 Hz o CA 220 V~ 50–60 Hz	Depending on version; see rear label
Consumption	Max. 400 VA	

Dimensions and weight

Item	Specification	Notes
Dimensions	330 (L) $\times$ 275 (W) $\times$ 589 (H) mm	
Net weight	18,1 kg	

Altitude

Item	Specification	Notes
Up to 2000 m	No purity reduction	
2000 m to 4000 m	Oxygen purity may be $< 90\%$	Varies according to air pressure within the altitude range of 0 m to 4000 m (fig. 2)

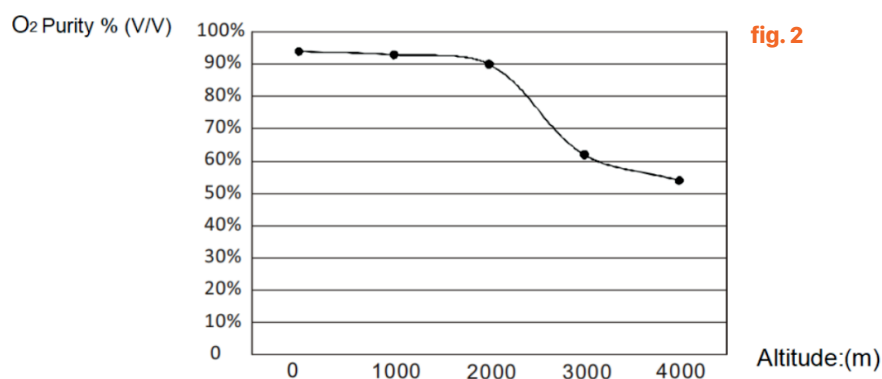Classification and Protection

Item	Specification	Notes
Class	Class II	At the compressor outlet
Applied part	Type BF	The applied part (cannula) is type BF
Operating mode	Continuous	
Degree of protection	IP21	

Storage and transportation conditions

Condition	Use (normal)	Storage/transport	Notes
Temperature	5 °C to 40 °C	-20 °C to +55 °C	
Relative humidity	20% to 80%	≤ 95% (non-condensing)	
Atmospheric pressure	86 kPa to 106 kPa	50 kPa to 106 kPa	
Position	Vertical	Vertical	The device must be stored, transported, and used only in the vertical position.

 **Caution:** if the device has been exposed to very high or very low temperatures, allow it to acclimate to room temperature (approx. 2 hours) before installation.

Information for technical service

- The manufacturer can provide circuit diagrams, component lists, descriptions, and calibration instructions to assist qualified service personnel in repairs.
- PCB fuse technical data:** T, 3,15 A, 250 V.

3. Unpacking and inspection


The oxygen concentrator is packaged to protect it during transportation and storage.


- Carefully remove the concentrator and all components/accessories from the box.
- Inspect the exterior surface for damage (impacts, dents, scratches, etc.).
- Inspect all components and accessories.
- If any damage is detected, contact your supplier.

Note: if the device will not be used immediately, keep the box and packaging materials.

4. Installation and operation

4.1 Schematic diagram and component functions

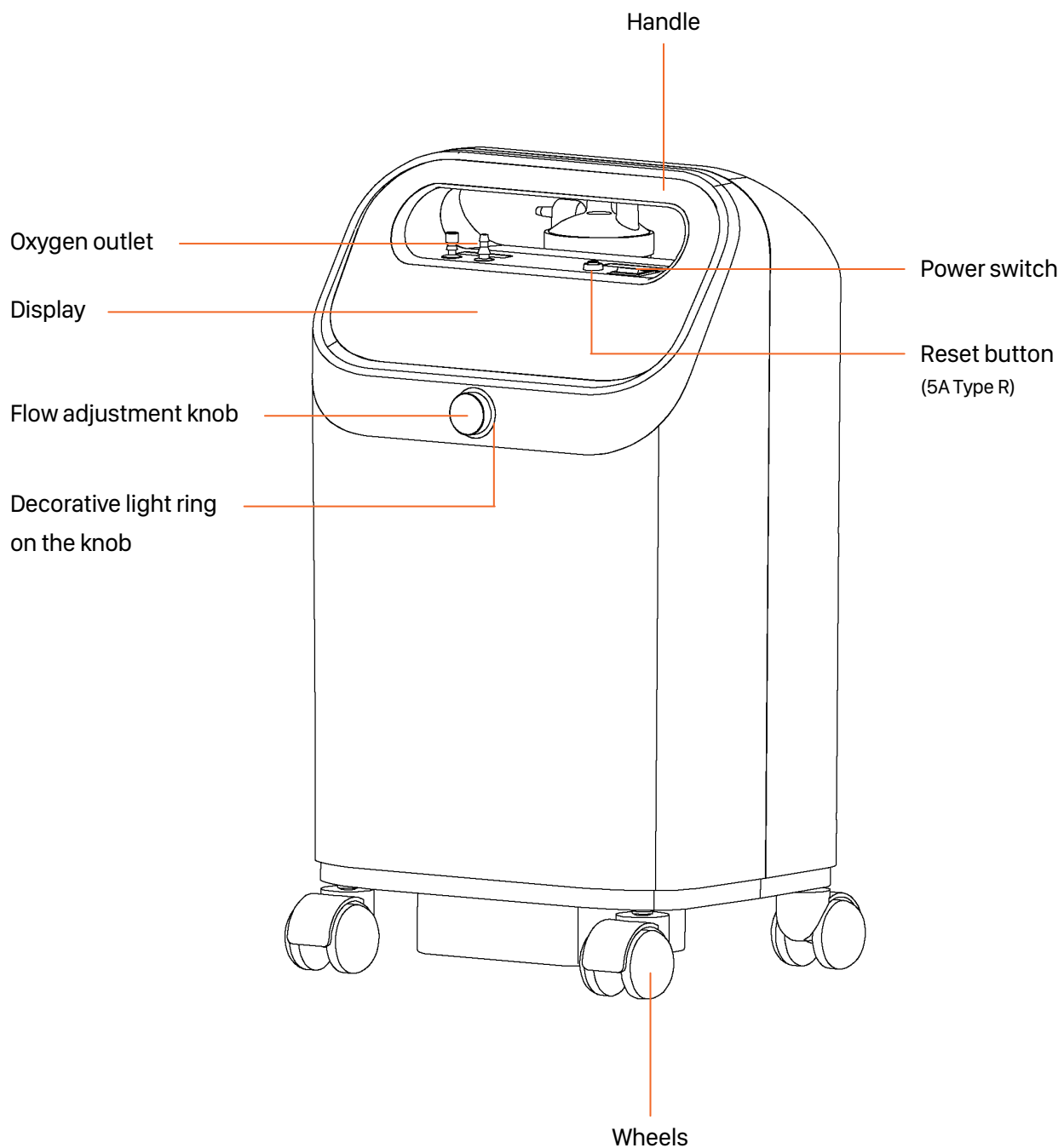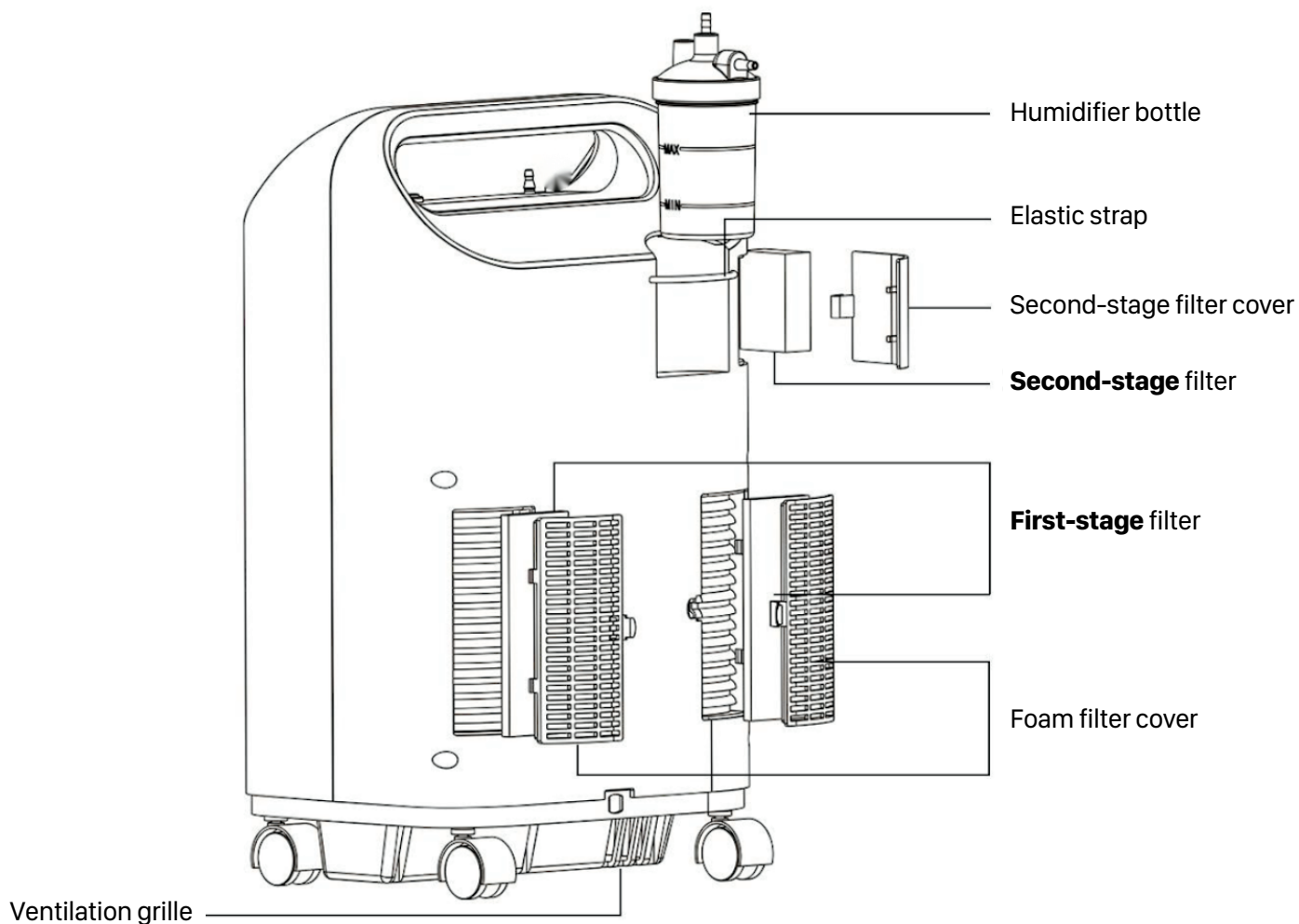Information for technical service

The reset button is a reusable circuit breaker. An ejected reset button (**fig. 3b**) indicates that the current in the circuit exceeds the rated current of the reset button; therefore, the oxygen concentrator power input is automatically disconnected and the device stops operating. For the normal state of the reset button, see **fig. 3a**

Note: if the reset button is ejected, first turn off the oxygen concentrator power switch, then unplug the power cord, and finally **contact a qualified supplier or the manufacturer.**

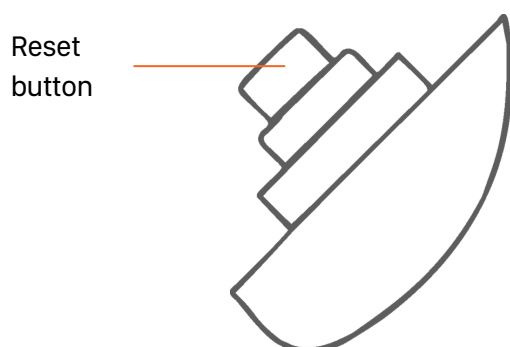
fig. 3a

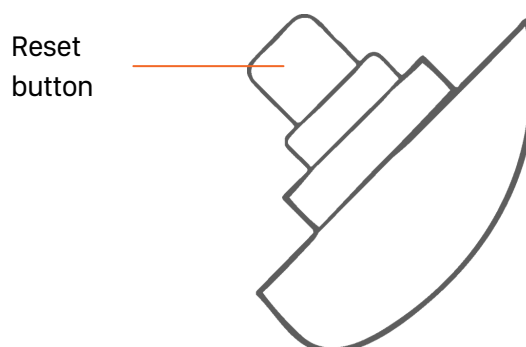
fig. 3b

4.2 Schematic diagram and use of the remote control

When the main power switch of the device is turned on, different functions can be activated by pressing the buttons on the remote control (see reference image).

1 Power ON/OFF

When the oxygen concentrator is in standby mode, press the power button  on the remote control to switch it to operating mode and start operation.

When the concentrator is operating, press the button again  to return to standby mode.

2 Increase Flow

Each time the **+** button is pressed, the oxygen flow increases by 0.1 L/min.

Flow adjustment range: 0.5 L/min to 5.0 L/min.

3 Disminuir flujo

Each time the **-** button is pressed, the oxygen flow decreases by 0.1 L/min.

Flow adjustment range: 0.5 L/min to 5.0 L/min.

4 Timer

Each press of this button increases the timer value on the display by **30 minutes**, and the operating time of the oxygen concentrator increases accordingly.

If the button is pressed repeatedly, the time changes cyclically in the following sequence:

00 h 30 min → 1 h 00 min → 1 h 30 min → ... → 8 h 00 min → continuous operation (timer off) → 00 h 30 min ...

If the button is pressed and held for **1.5 seconds**, the timer value increases rapidly following the same sequence. At the end, the device exits timer mode and returns to continuous operation mode (timer off), which is the normal operating mode.

5 Voice

Press the voice button once to activate the intelligent voice function. The device will announce **Voice on.**

Press the button again to deactivate the voice function. The device will announce **Voice off.**

6 Mute Alarm


When the oxygen concentrator is operating, press this button to mute all audible alarms. The mute symbol  will appear on the display, while all other information remains unchanged. Press the button again to reactivate the alarm sound and the symbol  will disappear.

 This function applies only to the current alarm state. If the alarm condition changes or is cleared, the button must be used again to manage the sound.



7 Help (SOS)

When the concentrator is operating, press this button and the device will continuously play the message **Help, please.**

Press the button again to stop the message. During SOS message playback, the concentrator continues operating normally.

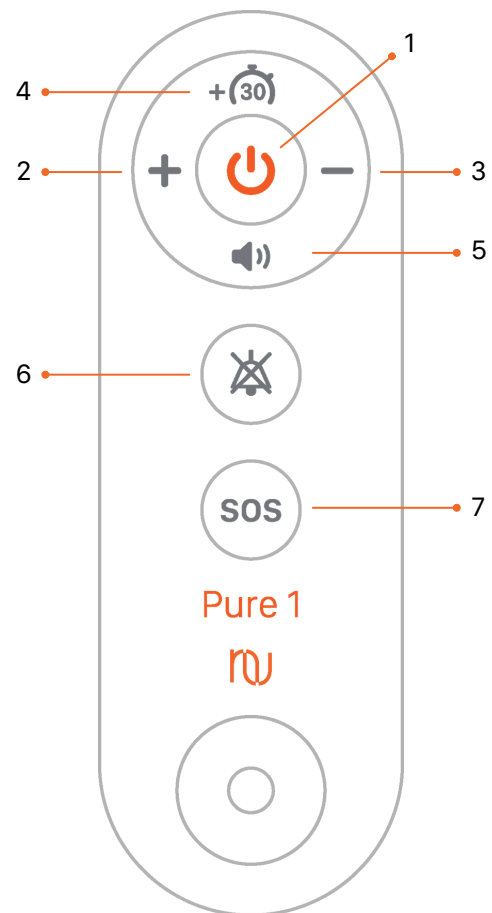4.3 Operation

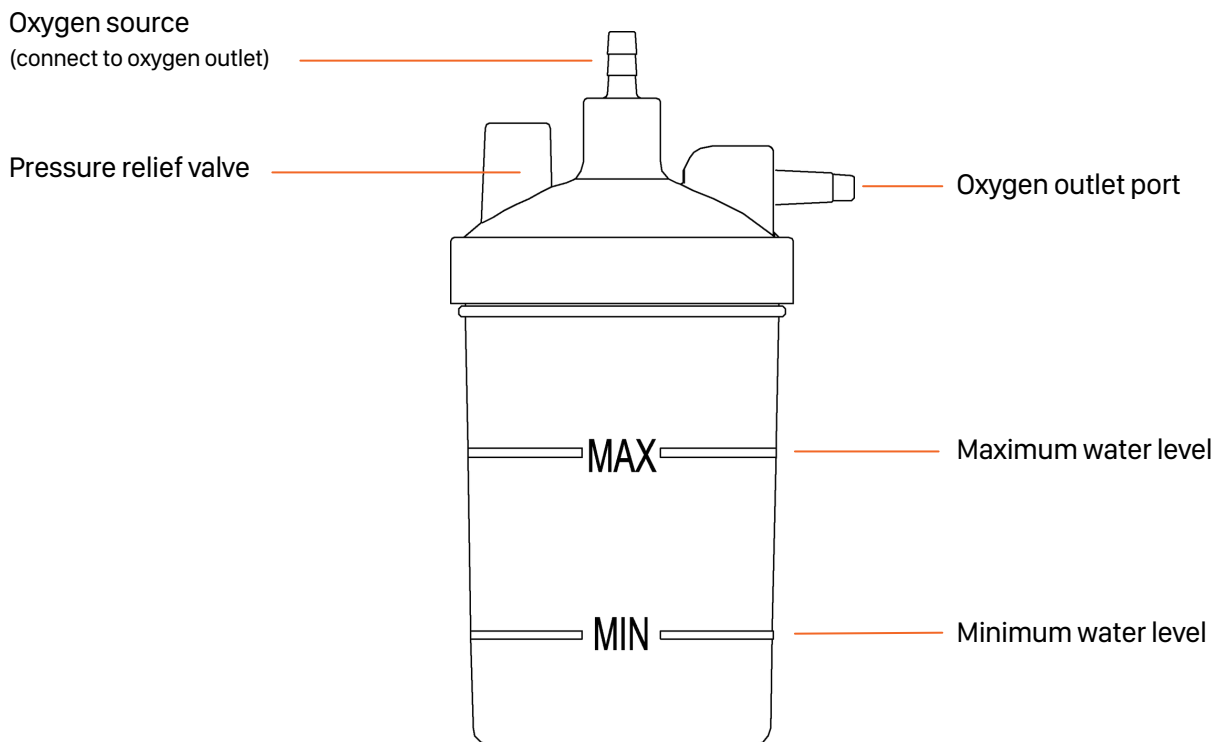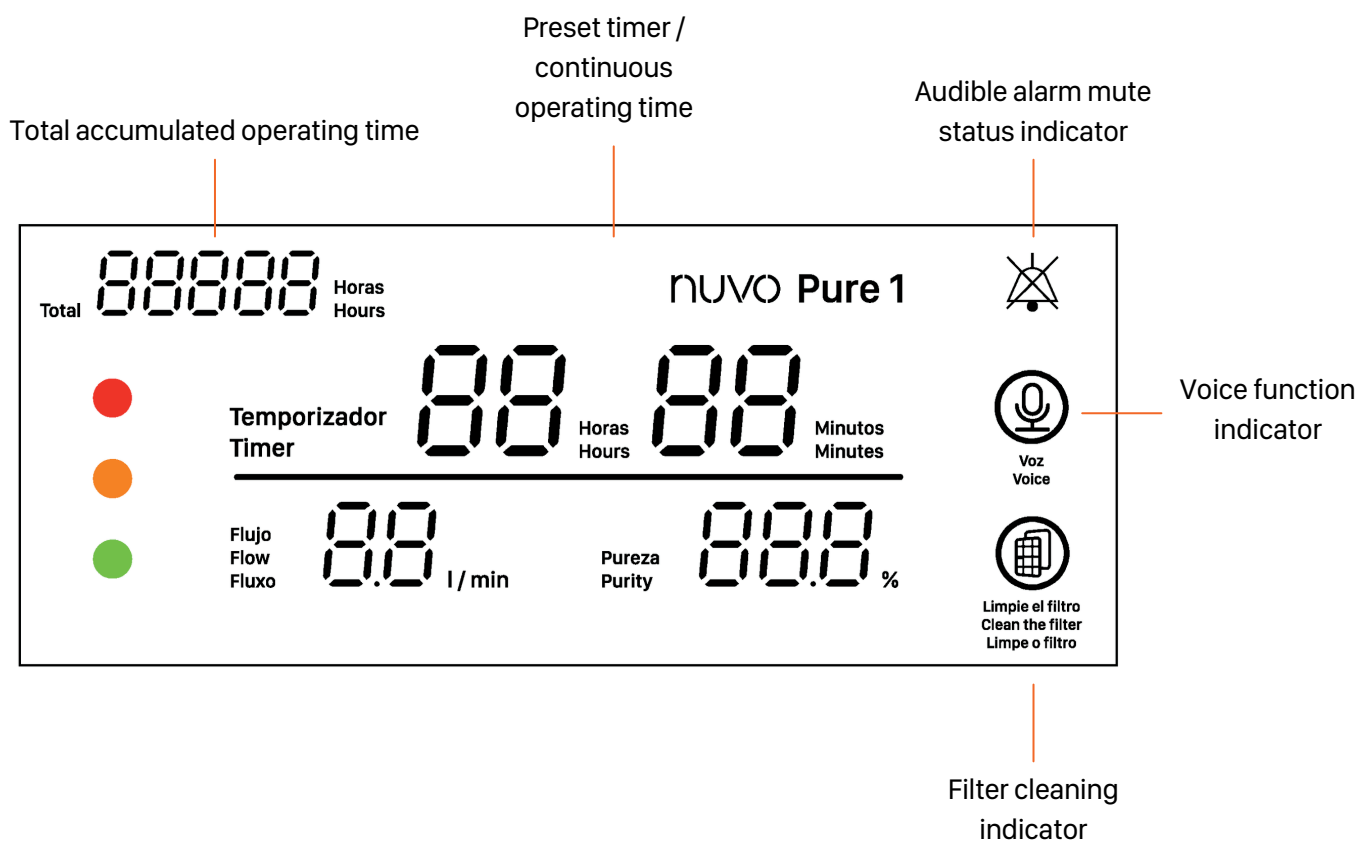Preparation

- 1 Place the concentrator on a stable, vertical, and well-ventilated surface.
- 2 Make sure the power switch is  (OFF).

WITH humidifier bottle (recommended if prescribed by a physician)

- 1 Remove the humidifier bottle upward (release the elastic band if applicable).
- 2 Unscrew and fill with purified or distilled water. The level must be between MIN and MAX.
- 3 Screw the cap back on tightly to prevent leaks.
- 4 Connect the tubing between the air inlet (oxygen source) of the bottle and the oxygen outlet of the concentrator.
- 5 Connect the silicone tube of the firebreak assembly to the outlet of the humidifier bottle.
- 6 Connect the nasal cannula to the other end of the firebreak assembly.

Important: replace the water in the humidifier bottle **before each treatment**.

WITHOUT humidifier bottle


- 1 Disconnect the nasal cannula and firebreak assembly from the humidifier bottle outlet.
- 2 Connect the nasal cannula directly to the oxygen outlet of the concentrator.

Power-on and flow adjustment


- 1 Set the power switch to  (ON).
- 2 After approximately 2 seconds, the display turns on and the green status indicator lights up. The following information is shown on the display: **flow, purity, individual operating time, and total accumulated time**.

Normal sounds:


- During operation, a “hiss... poof...” sound may be heard every few seconds (normal switching and exhaust).
- At power-on, a “smack” sound may be heard due to rapid internal valve switching (normal).
- Oxygen purity increases over time and reaches a stable value within 30 minutes.

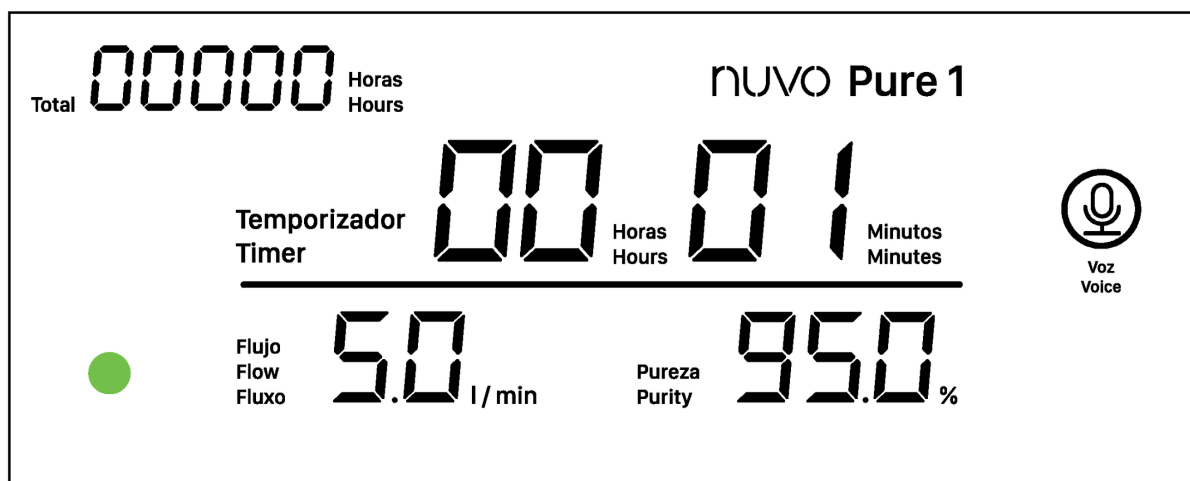- 3 Adjust the flow by turning the flow adjustment knob to the prescribed value. **Range: 0.5–5.0 L/min.** The flow rate and purity are displayed on the screen in real time.

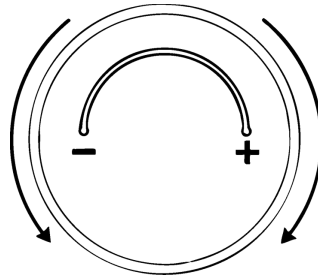- 4 Place the nasal cannula by inserting the tips into the nostrils, routing the tubing over the ears, and adjusting it under the chin.

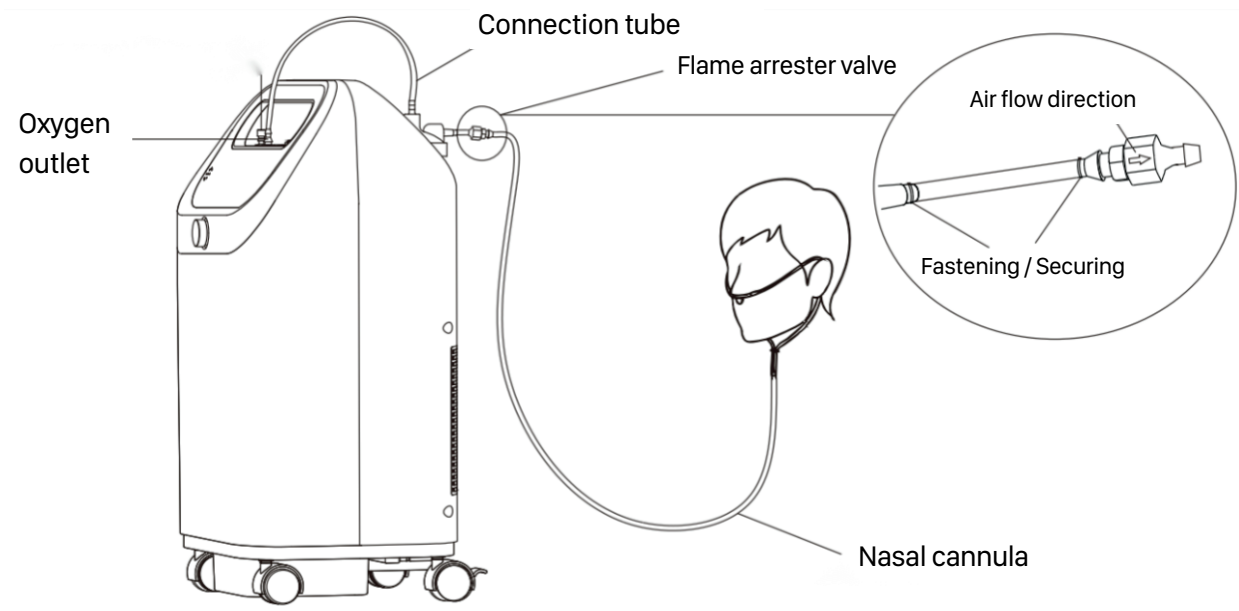Important:

- The required oxygen concentration is usually achieved within **2 minutes** after starting the device or adjusting the flow to a new value.
- To check oxygen output: you may submerge the tip of the cannula in water and observe bubbles, or feel the airflow with your hand.
- If there is no oxygen output, check for leaks in the connections or contact your supplier.
- **DO NOT** wear the cannula while adjusting the flow.
- **DO NOT** adjust the flow while receiving oxygen therapy. Adjust first, then place the cannula.
- The nasal cannula is **single-use only**. Dispose of it after use.
- **DO NOT** use if the package is **open or damaged**. Request a replacement from your supplier.
- Always follow medical instructions regarding inhalation time and flow rate.

Statement:

- Oxygen output settings must be adjusted individually for each patient and depend on device configuration and accessories. Proper cannula placement is essential for treatment effectiveness. Periodically reassess output settings to ensure continued effectiveness.

4.4 Connecting and using the mobile app

This oxygen concentrator can be paired with the **SmartLife** app to view device information and use certain functions from a mobile phone.

Downloading the app

Download the **SmartLife** app from the Apple App Store, Google Play Store, or the app store corresponding to your mobile device. You can also scan this QR code.

Registration or log in

Open the app and create an account using your mobile phone number or email address. If you already have a registered account, log in with your access details. The app may also allow quick or biometric login methods, depending on your phone's operating system.

- 1 Tap Log in with Password. Tap Register in the top right corner of the page to go to the registration page.
- 2 Enter your mobile phone number or email address, agree to the User Agreement, Privacy Policy, and Third-Party Information Sharing List, and then tap **Get Verification Code**.
- 3 Enter the verification code you received to go to the password setup page. Set a password as instructed and tap Done.


Requirements before pairing

Before starting the connection between the device and the app, check the following:


- 1 The concentrator must be **powered on and operating**.
- 2 The mobile phone must have **Wi-Fi and Bluetooth** enabled.
- 3 The app must have location permission, in addition to Wi-Fi and Bluetooth access, in order to correctly detect the device.

Pairing the concentrator with the app


- 1 On the main screen of the app, tap **Add Device or the (+) symbol**.
- 2 Keep the oxygen concentrator powered on and running.
- 3 On the concentrator interface, press and hold the Alarm Cancel/Mute button for **more than 10 seconds**.

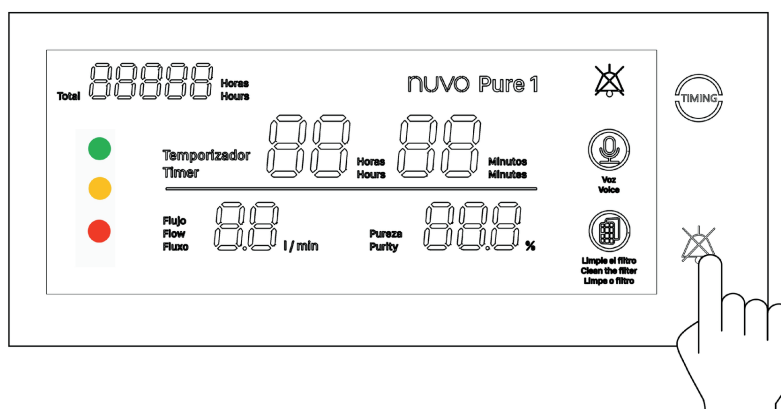- 4 When the device emits a single **"beep"** sound, it means it has entered pairing standby mode.
- 5 The app will automatically search for nearby devices.

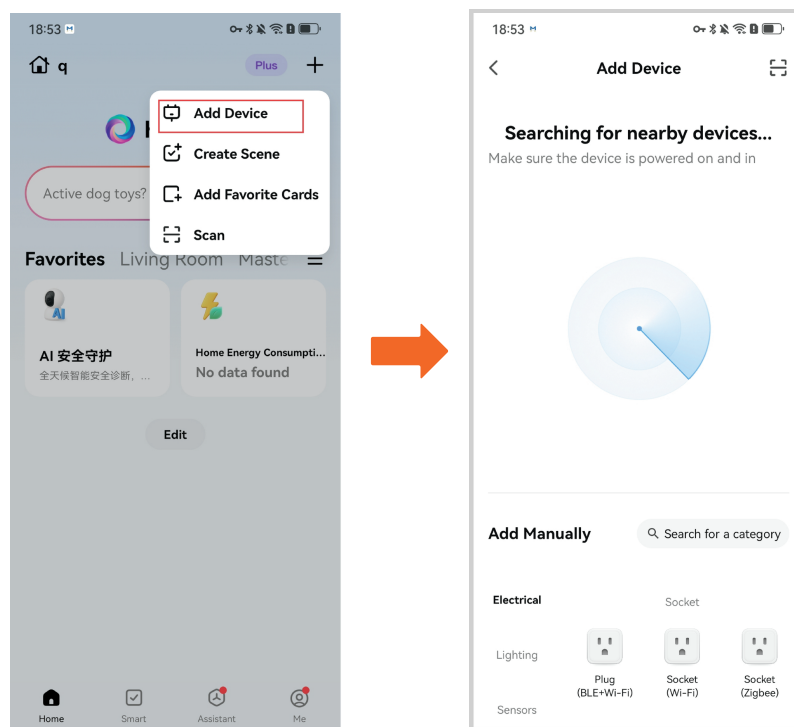- 6 When the device appears on the screen, select it.
- 7 Enter the Wi-Fi network information when prompted by the app and continue the process.

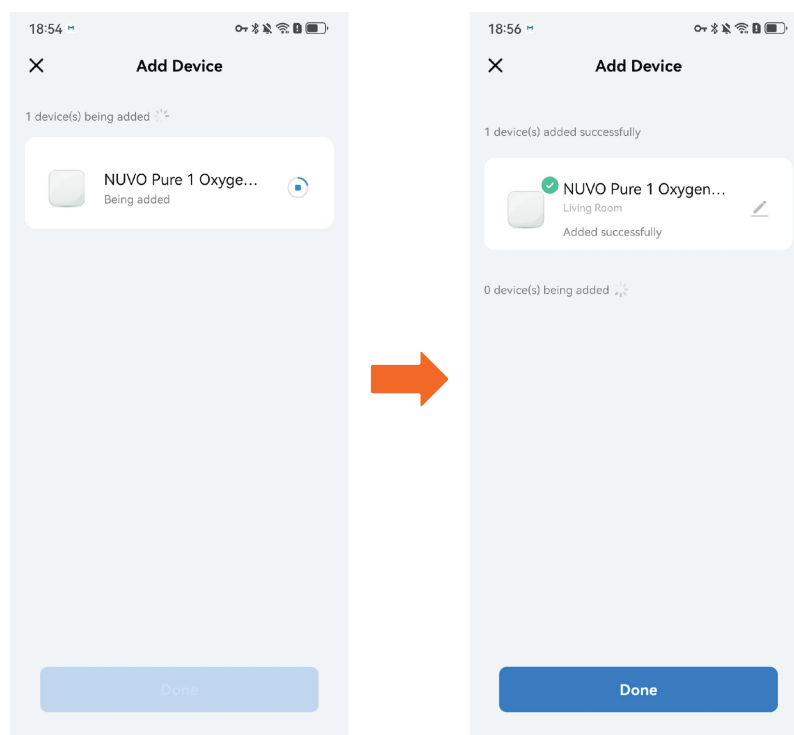- 8 Wait until pairing is completed successfully. Once the device is added successfully, tap **Done** to finish.

Functions available in the app

Once the device is paired, the app allows you to view and manage functions such as:

Notes:

- Use the app only as a support tool for operating the device.
- Always follow medical instructions regarding usage time and the prescribed flow setting.
- Do not change the flow adjustment without medical instruction. The flow must be set according to the corresponding prescription.
- If the device does not appear in the app, check that the concentrator is operating, that the phone has Wi-Fi, Bluetooth, and location enabled, and repeat the pairing procedure.
- This equipment is not entitled to protection against harmful interference and must not cause interference to duly authorized systems. For more information, please visit the Anatel website: www.anatel.gov.br

4.5 Alarms and indicators

This chapter describes the alarms and how to respond.

Power failure alarm test

With the device operating, unplug the power cord from the outlet. An audible alarm ("**beep beep beep, beep beep, beep beep beep, beep beep**") will sound, indicating that the alarm system is functioning properly.

Alarm priorities and indicators

Element	Low alarm	High alarm
Meaning	Requires caution and monitoring	Requires immediate response
LED	Solid yellow	Flashing red
Sound	"beep... beep..." (approx. every 20 seconds)	Repetitive pattern with groups of sounds and pauses
Display	Error code displayed in the center	Error code displayed in the center

In accordance with the requirements of the relevant standards, the volume of the audible alarm signal meets the specified requirements. The measured sound pressure range of the audible alarm signal is described below:

Alarm condition	A-weighted sound pressure level, average over the measurement surface (dB)	Remarks
High priority	70.59	Background noise: 31.9 dBA
Low priority	61.22	Background noise: 31.9 dBA

Muting the alarm

When an alarm sounds, press the  **mute/cancel** button to turn off the sound.
Press it again to resume the alarm sound.

Alarm code table (causes and responses)

Code	Priority	Condition	Device response	User action
E1	Low	Flow < 0.3 L/min for 30 s	Continues operating	Check connections, leaks, and settings. Ensure tubing/cannula is not blocked.
E2	Low	After 1 min 55 s, O ₂ purity ≤ 82% for > 2 s	Continues operating	Check filters, ventilation, and connections. Contact supplier if persistent.
E5	Low	Low mains voltage (out of spec)	Continues operating	Check the power supply. If the voltage is unstable, use a stabilizer/UPS as recommended by the supplier.
E6	High	Power failure	Stops operating	Restore power safely. If continuity is required, use a backup power source (UPS/battery) as recommended.

5. Cleaning and maintenance

5.1 Cleaning procedure

Before cleaning: turn off the device and disconnect it from the power supply.

Cleaning

- 1 Immerse the chamber components in a mild, pH-neutral soap solution for 5 minutes or more.
- 2 While immersed, brush channels and corners with a soft-bristle brush for at least 5 passes.
- 3 Remove and rinse with purified water for 2 minutes or more, ensuring all detergent residue is removed.

Drying

- 1 Allow to air dry naturally or dry with a sterile, low-lint cloth.

Note: if the device is not used for a long period, it is recommended to clean it using the above method before storage.

5.2 Maintenance

No special maintenance is required by the patient. Periodic maintenance is performed by the supplier to ensure reliable operation.

Preventive maintenance

- **First-stage filter:** clean every 200 hours or once per month.
- **Second-stage (HEPA) filter:** replace every 1000 hours (or more frequently depending on environment).

Oxygen concentration verification

- Verify concentration every 15,000 hours or 5 years. Must be performed by manufacturer-trained/certified technical personnel.
- Use original spare parts if replacement is required.

Alarm system verification (monthly)

- With the device operating, block the oxygen outlet. After approximately 40 seconds, the yellow light should turn on and the alarm should sound.
- Press the Mute audible alarm button to silence the sound; the bell icon appears on the display.
- Press again to reactivate the sound; the icon disappears.

Note: if failures occur that cannot be resolved, contact the supplier.

6. Troubleshooting

Problem	When it occurs	Possible causes	What to do
The device is operating, the yellow light is on, and the display shows "E1".	Flow is less than 0.3 L/min for 30 seconds.	- The first-stage and/or second-stage filter is blocked. - The oxygen outlet, tube, or humidifier bottle is blocked.	- Check the first-stage and/or second-stage filter. - If they are blocked or have been used for more than 1000 hours, replace the second-stage filter and wash the first-stage filter. ⚠ Do not use the device if the first-stage filter is wet or if the second-stage filter is not installed. - Check the oxygen outlet, tube, and humidifier bottle. Remove any obstruction.
The device is operating, the yellow light is on, and the display shows "E2".	After 1 min 55 s of operation, oxygen purity is $\leq 82\%$ for more than 2 seconds.	- The first-stage and/or second-stage filter is blocked. - The ventilation grilles are blocked.	- Check the first-stage and/or second-stage filter. - If they are blocked or exceed 1000 hours of use, replace the second-stage filter and wash the first-stage filter. ⚠ Do not use the device if the first-stage filter is wet or if the second-stage filter is not installed. - Ensure the ventilation grilles are not blocked and allow proper heat dissipation.
The device is operating, the yellow light is on, and the display shows "E5".	Mains voltage is 176 V or 184 V for more than 1 second.	Unstable or low voltage.	Use a voltage regulator.
The device stops, the red light flashes, and the display shows "E6".	Power supply is interrupted during operation.	- Power outage. - Circuit current exceeds the capacity of the reset button.	- Make sure the plug is securely connected. - Connect to an active power outlet. - Press the reset button. ⚠ Turn off the power before pressing the reset button. This operation must be performed by qualified technical personnel.
The device is operating and flow is properly set, but little or no oxygen is delivered.		Oxygen leak between the humidifier bottle and the cap.	Tighten the cap securely until there are no leaks.
		The pressure relief valve of the humidifier bottle is open.	Gently shake the bottle to close the valve.
		Incorrect connection of the oxygen tube.	Reconnect the tube correctly.
		Defective oxygen tube.	Replace the oxygen tube.

Note: if any of the problems persist or if another fault is detected that cannot be identified, contact the supplier.

Common user issues and solutions

Problem	Possible causes	What to do
Dry mouth, nose, or throat.	Using the device with the mouth open causes air to escape through the mouth and dry the airways.	Keep the mouth closed during use or use a full-face mask. Consult your physician for guidance.
Dry, cold, congested nose; runny nose; cold symptoms.	Cold airflow may irritate the nasal mucosa, causing dryness and inflammation.	Use the device with the humidifier bottle filled with water. Consult your physician and continue treatment unless advised otherwise.
Itching or redness in the nose.	Possible allergy to nasal cannula material.	Consult your physician and the equipment supplier. Use a nasal cannula made of non-PVC material.

7. EMC information (electromagnetic compatibility)

This oxygen concentrator complies with the EMC (Electromagnetic Compatibility) standard and **IEC 60601-1-2:2014+A1:2020** for the safe use of medical electrical equipment. The EMC standard is intended to ensure the safe use of medical electrical equipment. It specifies that interference caused by electromagnetic waves generated by other equipment, as well as electromagnetic interference emitted by devices such as mobile phones or other electronic equipment, must be kept within a specified range. **IEC 60601-1-2:2014+A1:2020** indicates that detailed information about the electromagnetic environment must be provided to ensure the safe operation of the equipment.

The EMC-related technical instructions are described below (for more details, refer to **IEC 60601-1-2:2014+A1:2020**). When this oxygen concentrator operates within the electromagnetic environment specified in these EMC technical data, its basic performance is not affected.

7.1 Definition of EMC

EMC (Electromagnetic Compatibility) is the ability to:

- 1 Not emit electromagnetic interference beyond the permitted limits that could affect nearby equipment (emissions).
- 2 Operate normally in an electromagnetic environment where other equipment generates interference (immunity).

7.2 EMC technical instructions

TABLE 1: Electromagnetic immunity · Manufacturer's guidance and declaration

This device is intended for use in the electromagnetic environment specified below. The user must ensure that it is used in such an environment.

Emissions test	Compliance	Electromagnetic environment · Guidance
RF emissions CISPR 11	Group 1	The device uses RF energy only for its internal operation. Therefore, its RF emissions are very low and are not likely to cause interference with nearby electronic equipment.
RF emissions CISPR 11	Class B	The device is suitable for use in all establishments, including domestic establishments and those directly connected to the public low-voltage power supply network that supplies buildings used for residential purposes.
Harmonic emissions IEC 61000-3-2	Class A	
Voltage fluctuations / flicker emissions IEC 61000-3-3	Complies	

- Failure to follow these guidelines may increase emissions or reduce the immunity of the equipment.
- Portable RF communication equipment (including antenna cables and external antennas) shall be used no closer than 30 cm to any part of the concentrator (including cables). Otherwise, performance may be degraded.
- If continuous operation is required during mains power interruptions, it is recommended to power the equipment using a UPS (uninterruptible power supply) or battery.

Nota: the power cord is NOT detachable; only service personnel may replace it in case of failure. Circuit diagrams may be provided by the supplier/manufacturer for maintenance by qualified technical personnel.

TABLE 2: Guidance and declaration · Electromagnetic immunity

This device is intended for use in the electromagnetic environment specified below. The user must ensure that it is used in such an environment.

Immunity test	IEC 60601 test level	Compliance level	Electromagnetic environment · Guidance
Electrostatic discharge (ESD) IEC 61000-4-2	±8 kV contact ±2 kV, ±4 kV, ±8 kV, ±15 kV air	±8 kV contact ±2 kV, ±4 kV, ±8 kV, ±15 kV air	Floors should be wood, concrete, or ceramic tile. If floors are covered with synthetic material, the relative humidity should be at least 30%.
Electrical fast transients / bursts IEC 61000-4-4	±2 kV for power supply lines ±1 kV for input/output lines	±2 kV for power supply lines	Power supply quality should be that of a typical commercial or hospital environment.
Surge IEC 61000-4-5	±0.5 kV, ±1 kV line to line ±0.5 kV, ±1 kV, ±2 kV line to earth	±0.5 kV, ±1 kV line to line	Power supply quality should be that of a typical commercial or hospital environment.
Voltage dips, short interruptions and voltage variations on power supply input lines IEC 61000-4-11	<5% Ut (>95% dip in Ut) for 0.5 cycle <5% Ut (>95% dip in Ut) for 1 cycle 70% Ut (30% dip in Ut) for 25/30 cycles <5% Ut (>95% dip in Ut) for 5 s	<5% Ut (>95% dip in Ut) for 0.5 cycle <5% Ut (>95% dip in Ut) for 1 cycle 70% Ut (30% dip in Ut) for 25/30 cycles <5% Ut (>95% dip in Ut) for 5 s	Power supply quality should be that of a typical commercial or hospital environment. If continuous operation during power interruptions is required, it is recommended to power the equipment from a UPS or battery.
Power frequency (50/60 Hz) magnetic field IEC 61000-4-8	30 A/m	30 A/m	Power frequency magnetic fields should be at levels characteristic of a typical location in a commercial or hospital environment.
Conducted RF IEC 61000-4-6	3 Vrms (150 kHz to 80 MHz) 6 Vrms in ISM bands	3 Vrms (150 kHz to 80 MHz) 6 Vrms in ISM bands	Portable and mobile RF communications equipment (including antenna cables and external antennas) should not be used closer than 30 cm (12 inches) to any part of the device, including cables specified by the manufacturer. Otherwise, performance may be degraded.
Radiated RF IEC 61000-4-3	3 V/m and 10 V/m, 80 MHz to 2.7 GHz. Test at 385 MHz to 5785 MHz. Specifications for enclosure port immunity to RF wireless communications equipment.	3 V/m and 10 V/m, 80 MHz to 2.7 GHz. Test at 385 MHz to 5785 MHz. Specifications for enclosure port immunity to RF wireless communications equipment.	If higher immunity test levels than those specified in Table 9 are used, the minimum separation distance may be reduced. Minimum separation distances for higher immunity test levels shall be calculated using the following equation: $E = \frac{6}{d} \times \sqrt{P}$ where: • P is the maximum power in watts (W) • d is the minimum separation distance in meters (m) • E is the immunity test level in V/m

Note: Ut is the AC mains voltage prior to application of the test level.

8. Other notes

8.1 Packing list

- **Oxygen concentrator:** 1 unit
- **Disposable nasal cannula:** 2 units
- **Humidifier bottle:** 1 unit
- **Connection tube:** 1 unit
- **First-stage filter:** 2 units
- **Second-stage filter:** 1 unit
- **User manual:** 1 unit
- **Quality certificate:** 1 unit
- **Firebreak assembly (with manual):** 1 set

Note: accessories must be oxygen compatible, designed for oxygen therapy, biocompatible, and comply with applicable regulatory requirements.

8.2 Accessories

- **NV001_NebulizationKit** – Nebulization kit
- **NV001_NasalCannula** – Nasal cannula
- **NV001_FireDamperValve** – Flame arrestor valve
- **NV001_HumidifierBottle** – Humidifier bottle
- **NV001_SiliconeTube** – Silicone connecting tube
- **NV001_PrefilterBlack** – Black filter – 1st stage (pre-filter)
- **NV001_HEPAFilter** – White HEPA filter – 2nd stage
- **NV001_RemoteControl** – Remote control

- Use only cannulas and accessories manufactured or recommended by the supplier. The use of improper accessories may affect performance, reduce treatment effectiveness, and void the manufacturer's liability.
- Do not pile long tubing at the head of the bed: it may wrap around the head or neck during sleep.
- Do not connect additional equipment unless recommended by the supplier or a healthcare professional.
- If parts are missing, contact the supplier.

8.3 Disposal (end-of-life)

Waste disposal method

All equipment waste (cannulas, filters, etc.) must be disposed of in accordance with local civil authority regulations.

Equipment disposal

The manufacturer considers environmental criteria, and a large proportion of the components are recyclable. Follow local ordinances and recycling programs. At the end of its service life, dispose of the equipment and packaging in accordance with local laws and regulations.

Materials and components that may contact the patient/operator

- **Housing/handles:** Flame-retardant ABS (PA-765(+))
- **Labels:** White PP film, ink
- **Flow control knob:** ABS plastic with aluminum coating
- **Silicone connection tube:** Silicone rubber compound
- **Oxygen outlet:** Stainless steel 304
- **Main switch/power cable:** PVC
- **Main power supply:** PVC
- **Firebreak assembly:** Flame-retardant ABS (PA-765(+)) and silicone rubber compound
- **Nasal cannula:** PVC
- **Humidifier bottle:** ABS/PP
- **First-stage filter:** Polyester
- **Second-stage/HEPA filter:** Polyester

Notice: the manufacturer reserves the right to change product technology and appearance without prior notice.

9. Limited warranty

Grupo Thera Inc. warrants that the device will be free from defects in workmanship and materials and will perform according to specifications for:

1 year for the main device / 3 months for all accessories

The warranty applies from the date of sale by the manufacturer to the distributor. If the product does not meet specifications, the manufacturer will repair or replace (at its discretion) the defective material or part.

This warranty does not cover damage caused by accident, misuse, abuse, alteration, or defects not related to materials or manufacturing. The manufacturer assumes no liability for economic losses, loss of profits, overhead costs, or consequential damages arising from the sale or use of the product.

Some jurisdictions do not allow exclusion or limitation of incidental or consequential damages; in such cases, the limitation may not apply. To exercise warranty rights, contact local authorized distributors or the manufacturer.

Commercial / Distribution Contact

Nuvo® is a brand of GRUPO THERA INC. Imported by: GRUPO THERA INC., 43rd Fl., Oceania Business Plaza, Tower 2000, Punta Pacifica, Panama City, Panama.

Brazil: Theramart Brasil Ltda – CNPJ: 45.784.439/0001-90 – Rua César Augusto Dalsóquio, No. 4255, Warehouse A, Room 30, Salseiros District, Itajaí – SC, ZIP Code 88311-500, Brazil.

Chile: Theramart SpA – Av. Presidente Kennedy 5600, Room 507, Vitacura, Chile.

Contact:

BRAZIL – parceiro@grupo-thera.com

Technical Responsible: Priscila de Gouveia Zamboni Lins – CRF-SC 20958

LATAM – info@grupo-thera.com

ANVISA Registration:

82692699001

Manual version:

2.0 / JUN 2026

The batch number, date of manufacture, and serial number are indicated on the product label.

 This device must not be modified without authorization from the legal manufacturer.